# INSTRUCTION MANUAL

Model SS-650

*Snore*  *Stopper*™

HIVOX® 

## How It Works

## Attention!

- i. The background noise should be below 55 db while you use Snore•Stopper.
- ii. Normally, you might need a few days to get used to this device, after this period you will start to enjoy it.

### P.S.

What you get:

One set of Snore•Stopper includes:

- \* Main unit of Snore•Stopper(1)
- \* Wrist-belt with magic tape and hook (1)
- \* 1.5 V AAA battery (1)
- \* Pack of self - adhesive electrode pads, 3 pair (1)
- \* Manual (1)
- \* Warranty card (1)
- \* Storage case (1)
- \* Screwdriver (1)

## Contents

1. Introduction .....	2
2. Names of the components .....	3
3. Operation the unit .....	4
4. Maintenance and Trouble-shootin.....	8
5. User FAQs .....	10
6. Caution .....	13
7. Warning .....	14
8. Technical Specification .....	15

## Introduction


At this point we thank you for buying Snore•Stopper. It was a good choice. But firstly please carefully read the manual to get used to the right operation and the possibilities of the equipment. This way you will be able to enjoy the benefits of Snore•Stopper and you will have pleasure with it for a long time.

Our modern life is tensed and stressful. Therefore one of the most important things is a silent night and the possibility to recuperate during sleeping. In estimation, about 30% to 40% of the middle aged of the current population snore, and more and more people are affected by snoring! Many people think snoring is only a little problem and it will not influence your health. This is a big mistake! It might cause high blood pressure, increase the risk of heart disease, stroke, fatigue, depression, etc. But one thing is for sure: even in minor severe cases of snoring, we annoy our beloved partners and keep them from recovering during the night!

But now all these can be stopped by using Snore•Stopper! Snore•Stopper is a new device, designed to reduce snoring by a biofeedback process: when the device detects snoring, a small electronic gentle pulse will be sent automatically to stimulate the snorer's skin through nerve for around 5 seconds until the snorer reduce the snore by the laryngeal muscles get back into a normal condition, or alter the sleeping position. This unique process results in a subconscious reaction, without waking the snorer up and without disturbing the brain that controls his sleep. You will see the great result when you go on using Snore•Stopper regularly for 6 to 8 weeks. Snore•Stopper is a safe and intelligent choice to reduce your annoying snoring, and lets your family expect a higher sleeping quality in the future!

Use only as directed. If any unusual symptoms appear, discontinue the treatment and see your doctor.

## The Components
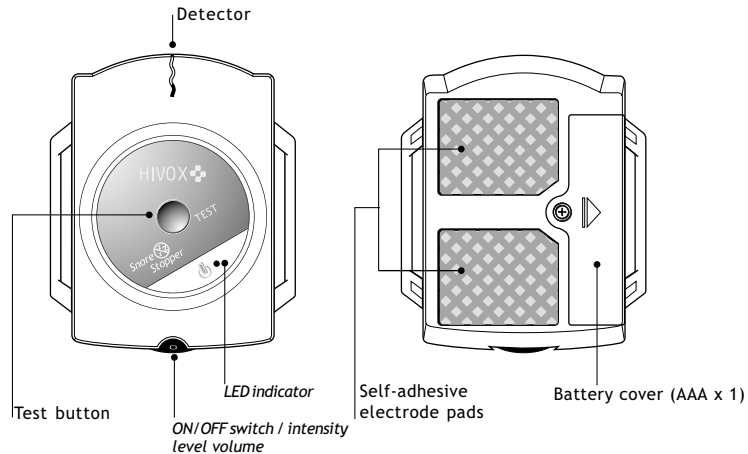
## Operate the unit

### I. Assembly before use

#### 1. Battery loading:


Please open the battery box with a screwdriver to load AAA battery. Make sure that the AAA + pole is placed towards "ON/OFF switch / intensity level volume" as shown in picture I.

2. **Then close the battery cover and fix it with the screw.**

### II. Stick the electrode pads

1. Open the vacuum pack and take out one pair of self-adhesive electrode pads. Then carefully detach the big, square plastic film on which the words "SELF-ADHESIVE ELECTRODE PADS" are printed (At this time, do not remove the other unprinted big plastic film.) as shown in picture II.
2. Stick the electrode pads, on the backside of the device, onto the black conductive rubbers. Gently press the pads with your finger onto the device, then take away the big plastic film as shown in Picture III.


(Tip: Whenever you remove the old electrode pads and put on new ones, press firmly with a thumb on the pads attaching them to the device, and move in a circular motion to ensure better adhesion.)
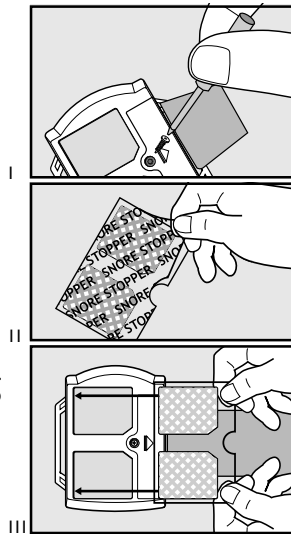
## Operate the unit

### III. Adjust the wrist-belt




1. After you have attached the electrode pads, you prepare to wear the Snore•Stopper. As it is shown in picture IV the wrist-belt has to be put into the left wing. Please fix the magic tape on this side now.
2. Now please remove the big plastic film and place Snore•Stopper 1 inch above the wrist as shown in picture V and hook the belt into the other "wing" of the device. Then tighten it and fix this end of the belt as shown in picture VI. (Do not pull the belt too tight; otherwise it will lead to poor circulation. In fact, the device will function well, when it touches the skin gently.)
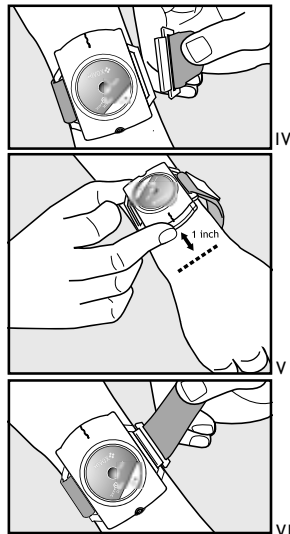
## Operate the unit

### IV. How to use the device

1. Test the function of the device before use. Switch on the power and then adjust the intensity to level 1. Check if the LED indicator (green light) is blinking.
2. Press the test button three times in a row. The device will release a weak electrical impulse that lasts for about 5 seconds. Now, you can adjust the intensity to the level you like.
3. Due to the variations of the skin sensitivity, you might not feel the electrical impulse. In that case, press the test button 3 times in repetition, and adjust the intensity to the level you like. By the way, the nerves are stimulated at nighttime easier. Therefore we advise, that you reduce the intensity for the night for 1-2 levels to have the best effect.
4. **Meaning of the LED light:**  
When you wear Snore•Stopper on your arm the LED normally changes from "red" to "green", which means, the device has contact to your skin properly and you can start using it. If the LED still shows red light, the device has no proper contact to your skin. Therefore, you should adjust the wrist belt into a position until you see that the LED indicator

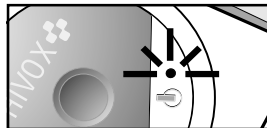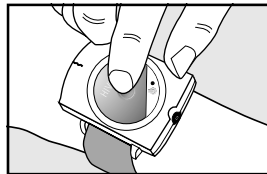
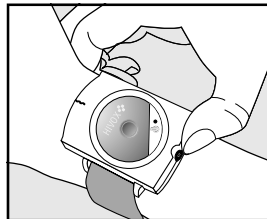
## Operate the unit

shows green light. If the device still does not work, please look at the MAINTENANCE & TROUBLESHOOTING section on page 8.

5. Our device's patented intelligent design can detect snores above 65 dB automatically, but without other background noise. By now, you can stop the annoying snores and let you and your family expect a better sleeping quality in the future!






### V. After Use

1. After use, please see the following instructions to prevent damage of the self-adhesive electrode pads: Unhook the wrist-belt and detach the device from your wrist gently. After you put off the device, please stick the big plastic film back onto the electrode pads. This can prevent the pads of contamination with dusty particles and will enhance their lifespan. This device will automatically switch off after 8 hours. In order to enhance the battery's life, please switch power off after you wake up.
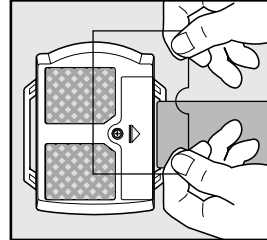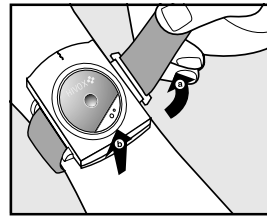
## Maintenance and Trouble-Shooting

1. The self-adhesive electrode pads should be kept clean. In order to maintain the device's performance, **DO NOT** rub the electrode pads with your hands, facial tissue, or wool.
2. The electrode pads tend to lose their adhesion over time, therefore weakening the intensity level. Moisten finger with a few drops of water and rub it gently against the electrode pads to improve contact. We recommend you replace the electrode pads after 30 days of usage (varies among individuals).
3. If you still do not feel the signal, please replace the old pads by a pair of new ones. Please refer to section 2. electrode pads sticking on page 4. (Note: Generally, the electrode pads can be used for 30 nights continuously.)
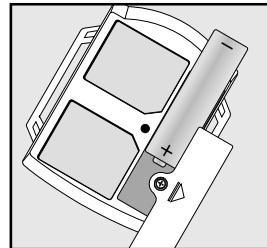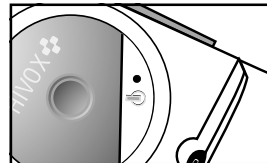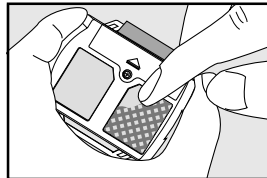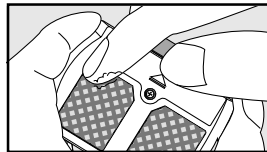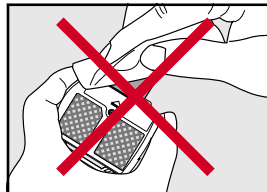
## Maintenance and Trouble-Shooting

4. If the LED indicator does not light during the testing procedure, please check the following:
  - \* Is the battery installed correctly? ("+" sign faces upward?)
  - \* Is the battery's energy low? If so, please replace by a new one.
5. If the points mentioned above are checked, and the device still does not work correctly, please send it back to your local dealer for repair. You are urgently reminded to include the warranty card!

**Note:** Always use the specified accessories in the manual. The use of other parts that are not approved by the manufacturer may cause hazard. Do not try to repair the device! There are no user serviceable parts inside.

## User FAQs

### 1. Q : Does Snore•Stopper really work ?

A : Yes, it works, but not for those with serious medical conditions as we indicate in our introduction and warning sections. Use only as directed. If any unusual symptoms appear, discontinue the treatment and see your doctor. The proper use of Snore•Stopper requires self-discipline, careful application and time to train reflexes.

We employ our special in-built electronic biosensor, which is designed to detect snoring, based on the biophysical properties of snores. Besides, due to its broad detection capability, no matter if your hands are under or outside the blankets, if the snores are above 65-70 dB and occur 3 times in a row, they will be detected and trigger the device to response automatically.

### 2. Does it really reduce my snoring?

A: Snoring is caused by a narrow airway. Firstly, when we sleep, our nerves are in a relaxed state, and when the air moves into the narrow airway during breathing, it causes the soft tissues of the throat (tonsils, soft palate & uvula) to vibrate. This vibration is the sound of snoring. Generally, snoring can be reduced if the snorers change their sleeping positions or are given a mild nerve stimulation. Based on this concept, Snore•Stopper is developed. When the snores are detected, the device will output a 5 seconds' impulse onto superficial skin, then snoring may be reduced in the following ways: i.) Snoring persons may change the sleeping position ii.) increase the muscle tense of the throat muscles due to nerves stimulation.

## User FAQs

### 3. What is dB ? How do I know whether Snore•Stopper can be used in my home?

A : dB (decibel) is an index of measuring the intensity of sound. If you are in an environment that has a sound intensity of above 70 dB continuously, we talk about a noisy background, normally, you will not live in this environment. The design of Snore•Stopper aims on the source of indoor background noise. You may read the following table to see which is the optimal sleeping condition for use.

Indoor background noise	Snore•Stopper Sensor
Above 65 dB (within 50 cm from TV)	Seldom Activated
60 dB (within 1 m from TV)	Seldom Activated
55 dB (within 1m from air-conditioner)	Seldom Activated
Below 55 dB	No Activated

Remarks : Under normal conditions, your home's indoor background noise will be 30-50 dB

### 4. Why do I feel the intensity so weak during the testing procedures ?

A: This may be due to three ways, i.) the difference in sex, age, weight and skin impedance will cause the variations of your feeling. In fact, the intensity of the output impulse lies within the safe range, so generally you should have the feeling. If your attention is called away, you may also have the low-intensity feeling. ii.) the skin's nerves are more sensitive during the nighttime. Since the device is designed to be used at night, you may have the low-intensity feeling when you test the device in daytime. iii) check the conductive electrode pads to see if you forget to attach them, or if they should be replaced due to extended use.

## User FAQs

### 5. When I 'simulate' the snores in a silent background, why does the device not work?

A: Since Snore•Stopper's special biosensor is designed according to the biophysical properties of snoring, so the device will react, and you may try to 'simulate' in such ways: say a long 1, 2, 3 or produce 3 consecutive snores with the same frequency and at least above 65 dB. Please try it again!

### 6. Why do I feel a special long time stimulation? Is something wrong with the device?




A: This is a special built-in design. As the microprocessor monitors your snoring frequency, if the device detects 15 times of outputs, it emits a 25-seconds long stimulation is to train the heavy snoring person to change the sleeping position.

### 7. Are the two sides of the electrode pad different? Does it matter if there is no proper contact?

A: Yes, the two sides are different. To tell them apart, the one side of the electrode pad, which must be attached to the device, has stronger adhesion (the side with printed words on the waxed protective film).

The device will work as long as the electrode pads are properly attached to the black conductive rubbers on the device, otherwise, electrode pads also can be use for last a while. Do not to place the two electrode pads too closely to each other, as this will affect conductivity. In the event the pads were stacked together, gently move them apart and resume proper attachment.

## Caution

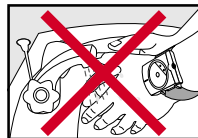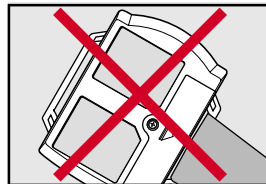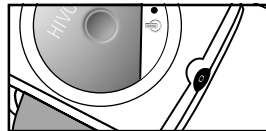1. The intensity level volume should be adjusted to ZERO before using.
2. DO NOT wear the device if the self-adhesive electrode pads are not stuck to it.  
Also, when you are pressing the test button, DO NOT touch the conductive rubber or electrode pads with your fingers.
3. DO NOT wear the device on the underside of the wrist.
4. DO NOT wash your hands when you are wearing the device. It is not waterproof.
5. When you are wearing the wrist-belt, DO NOT pull the belt strongly, otherwise the lifespan of device will be shortened.
6. DO NOT use with other medical electronic devices AT THE SAME TIME.
7. The device should be kept out of the reach of children.
8. Batteries must not be disposed in household waste. Return

them to public collection points or shops selling batteries of the same kind.

## WARNING !

1. This device is NOT suitable for persons with implanted cardiac pacemakers, and those monitored by electronic equipments (e.g. ECG monitors and ECG alarms)
2. Heart disease patients SHOULD NOT use this device except under advice of a medical doctor.
3. Snore●Stopper is not suitable for Sleep Apnea patients.
4. Do not use this device during pregnancy except with a medical doctor's consent.
5. Patients with acute diseases and infectious diseases may use this device only under a medical doctor's advice.
6. Patients with severe skin diseases or sensitive skin, or those in a state of bleeding SHOULD NOT use this device.
7. DO NOT use this device on over swollen, infected, or inflamed areas or skin eruption, e.g., phlebitis, thrombophlebitis, varicose veins, etc.
8. DO NOT place the device in the mouth, neck and eyelid region.
9. Operation in close proximity (e.g. 1 meter) to a short wave or microwave therapy equipment may produce instability in the device output.

## Technical Specification

Pulse amplitude:	Fixed 180µA peak into 1,000Ω load
Pulse rate:	Fixed 0.2 Hz
Pulse width:	Fixed 800 µs, Voltage adjustment 2.5-14V, peak into 1,000Ω load
Wave form:	Asymmetrical bi-phasic square pulse
Power source:	1.5V Volt AAA battery
Detective timer:	8 hours auto off
Operation environment:	10°C ~ +40°C, 30 ~ 85 % RH
Storage environment:	-20°C ~ +50°C, 10 ~ 95 % RH
Transport environment:	-10°C ~ +50°C, 35 ~ 85 % RH
Max wrist belt length:	22 cm, if you request a longer wrist belt, please contact your local dealer for replacement. ( All electrical specifications are ±20% 1,000Ω load )

*\* Design and specifications are subject to change without notice.*

**Note :** This device complies with the Electromagnetic Compatibility Requirement of EN 60601-1-2 as specified in EEC Directive 93/42/EEC.

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©2002 by HIVOX BIOTEK INC.